

Arterial Catheterization Sets

TRADITIONAL SELDINGER TECHNIQUE



SMOOTH NEEDLE PLACEMENT

- ▶ Low-friction needle with improved bevel requires less penetration force during insertion.
- ▶ Offers better control of needle placement within the vessel.



EXTENSION SET

- ▶ Integral extension set reduces potential clinician exposure to blood-borne pathogens.
- ▶ Provides less manipulation at the insertion site.



SUTURE WINGS

- ▶ The low-profile, angled suture wing improves patient comfort and supports the catheter insertion angle for optimal positioning.



IMPROVED GUIDEWIRE

- ▶ Depth markings to optimize positioning on spring-wire guides (straight wires only).
- ▶ Shorter length for improved handling.
- ▶ New: spring-wire guides longer than 35cm with convenient advancer.



HEMODYNAMIC MONITORING

ARROW
INTERNATIONAL

SMOOTH PLACEMENT


The polyurethane hydrophilic coated material facilitates smooth insertion, while the radius tip provides a smoother, less abrupt edge for an improved guidewire-to-catheter transition.

REDUCE COMPLICATIONS FROM KINKING AND VESSEL TRAUMA

- ▶ Arrow's coated polyurethane material is bio-compatible.
- ▶ Specifically designed to provide consistent catheter performance and resist waveform dampening, resulting in more reliable patient monitoring.

ADDITIONAL FEATURES

- ▶ ISO color-coded needle hubs enable easy recognition of the needle gauge size.
- ▶ Needle hubs are "V" notched to indicate bevel-up orientation.
- ▶ The versatile juncture hub is designed for multiple methods of securement.
- ▶ PVC-free; DEHP-free.
- ▶ New compact packaging of the spring-wire guide inside a catheter guard.
- ▶ Spring-wire guides, that are straight on both ends, have a reference mark that indicates when the tip of the spring-wire guide is at the needle tip.

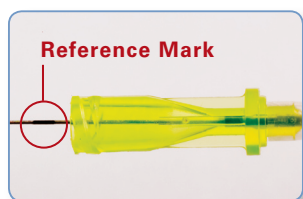

ARTERIAL CATHETERIZATION SETS

Product No.	Description	Needle Size	SWG Size	CE	sets/case
SAC-00818	18ga x 8cm	18ga x 5cm	.025" x 33cm - S	✓	10
SAC-01218	18ga x 12cm	18ga x 5cm	.025" x 33cm - S	✓	10
SAC-01618	18ga x 16cm	18ga x 7cm	.025" x 45cm - J	✓	10
SAC-02318	18ga x 23cm	18ga x 7cm	.025" x 60cm - J	✓	10
SAC-00520	20ga x 5cm	20ga x 4cm	.021" x 35cm - S	✓	10
SAC-00820	20ga x 8cm	20ga x 4cm	.021" x 35cm - S	✓	10
SAC-01220	20ga x 12cm	20ga x 7cm	.021" x 35cm - S	✓	10
SAC-01620	20ga x 16cm	20ga x 7cm	.021" x 50cm - S	✓	10
SAC-00522	22ga x 5cm	22ga x 4cm	.021" x 35cm - S	✓	10
SAC-00822	22ga x 8cm	22ga x 4cm	.021" x 35cm - S	✓	10
SAC-01222	22ga x 12cm	22ga x 4cm	.021" x 35cm - S	✓	10

Latex Free

"J" denotes spring-wire guide with J tip on one end and straight tip at the other end

"S" denotes spring-wire guide with straight tip on both ends

EACH SET CONTAINS:

- 1 Indwelling catheter with low-profile suture wing and extension tube with slide clamp for blood containment
- 1 Spring-wire guide - various sizes are available in straight and "J" tip (with straightener) configurations
- 1 Introducer needle

Caution: Rx only. Contents of unopened, undamaged package are sterile. Disposable. Refer to package insert for current warnings, indications, contraindications, precautions, and instructions for use.

To locate an Arrow Distributor or for ordering information, contact Arrow Customer Service at:

United States

Tel: 800.523.8446

Fax: 800.343.2935

International

Tel: +1.610.655.8522

Fax: +1.610.655.8566

www.arrowintl.com

ARROW
INTERNATIONAL

Distribution Worldwide:

Arrow offices are located in Belgium, Canada, Czech Republic, France, Germany, Greece, India, Italy, Japan, Mexico, Netherlands, Slovakia, South Africa, Spain, United Kingdom and the United States.

Arrow International, Inc., 2400 Bernville Road,
Reading, PA 19605 U.S.A. Tel: 610.378.0131

www.arrowintl.com

ENG-SAC-C 05/2007 CRB-07-0148

© 2007 Arrow International, Inc. U.S.A. All rights reserved.

